



Medtronic 社製 EnRhythmP1501DR 植込み後遠隔期に胸部違和感で受診

70 歳代男性。2010 年 7 月に間欠性房室ブロックに対して Medtronic 社製 EnRhythmP1501DR、心房リードは右心耳にメドトロニック社製 CapsureSP4592、右心室リードは中位中隔に日本ライフライン社製シンライン 438-35S を植込み。植込み後安定していたが、2012/5 に胸部違和感と脈の乱れを感じ受診。その際にとられた心電図。(補足：心房ペーシング、センシングはできていそう。心室リードが中隔のためペーシングでもⅡ誘導でも以前から narrowQRS、ペーシング正常。ただしいつもは自己房室伝導リズムだった) 原因と対処を議論したい。

点検時の心電図

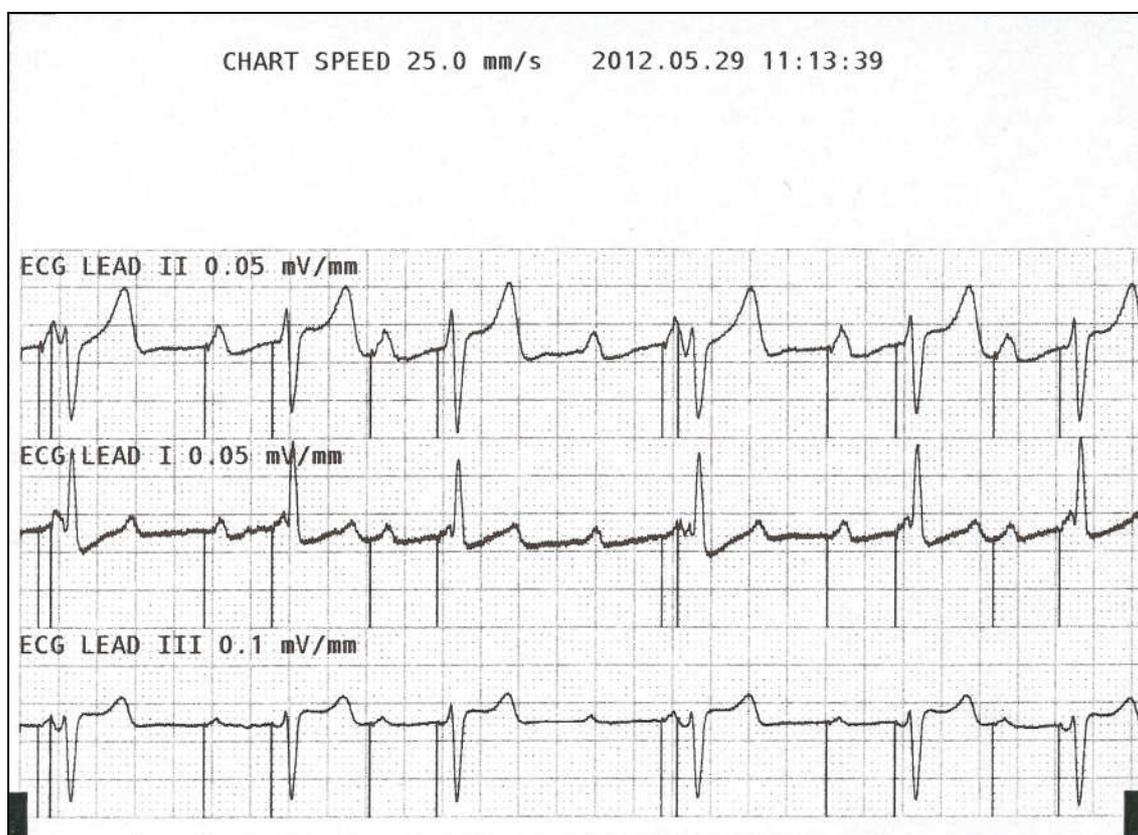

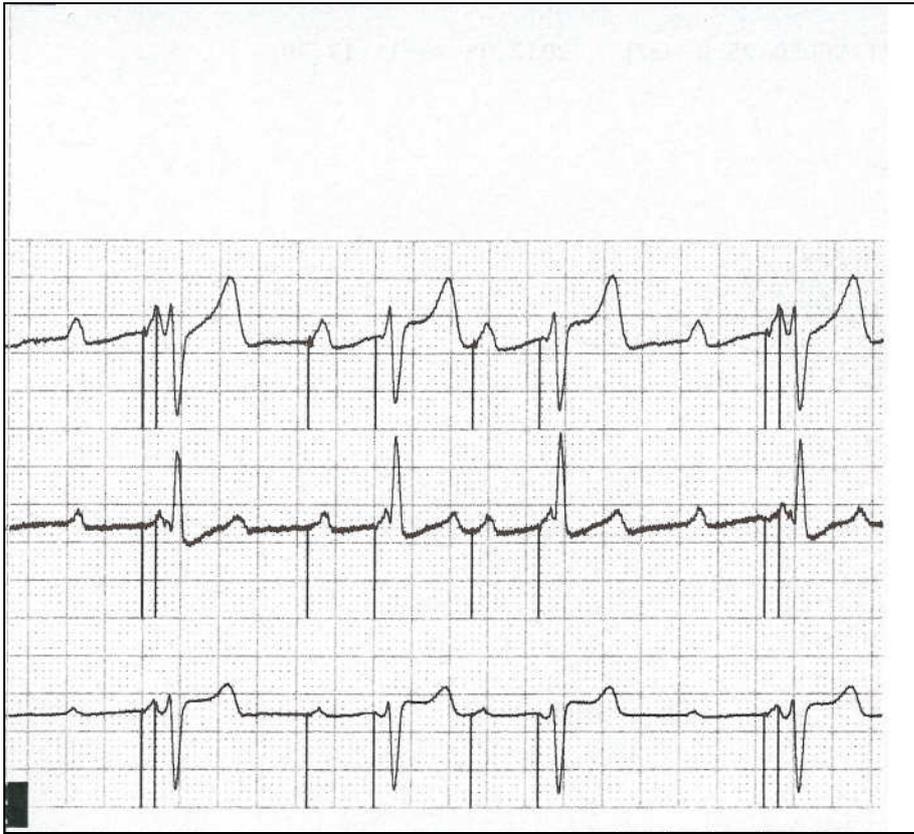

来院時のパラメータ

Pacing Summary			
Mode	Rates	AV Intervals	
Mode AAI<=>DDD	Lower 70 bpm	Paced AV 340 ms	
Mode Switch On	Upper Track 130 bpm	Sensed AV 310 ms	
	Upper Sensor 130 bpm		
Pacing Details			
Rate Response		Atrial	RV
Rate Response	7	Amplitude	2 V 2.5 V
Activity Threshold	Medium Low	Pulse Width	0.4 ms 0.4 ms
Activity Acceleration	30 sec	Sensitivity	0.3 mV 0.9 mV
Activity Deceleration	5 min	Pace Polarity	Bipolar Bipolar
		Sense Polarity	Bipolar Bipolar
AV Therapies		Arrhythmia Interventions	
Rate Adaptive AV	Off	A. Rate Stabilization	Off
Refractory/Blanking		A. Preference Pacing	Off
PVARP	250 ms	Post Mode Switch	Off
PVAB Interval	150 ms	V. Rate Stabilization	Off
PVAB Method	Partial		
A. Blank Post AP	200 ms		
A. Blank Post AS	100 ms		
V. Blank Post VP	200 ms		
V. Blank Post VS	120 ms		
Additional Features			
Non-Comp Atrial Pacing	On		
NCAP Interval	300 ms		
PMT Intervention	Off		
PVC Response	On		
V. Safety Pacing	On		